


Dried graminacee

Milestone Application Note for Acid Digestion MAXI-Food-056

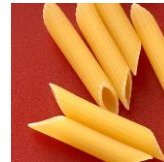

Summary

This method represents a guideline for the microwave acid digestion of Dried graminacee samples with the Milestone ETHOS UP and ETHOS EASY systems. This procedure may need to be optimized to match your specific analytical goals.

Instrumentation

Milestone ETHOS UP or ETHOS EASY microwave units, MAXI-44 high-throughput rotor, T2 infrared temperature control. Option P2 pressure control.



Sample weight and reagents

Sample weight	Up to 200 mg
Reagents	7 mL HNO ₃ 1 mL H ₂ O ₂

Microwave program

Step	Time	T2	Power
1	00:20:00	200°C	1800 W*
2	00:15:00	200°C	1800 W*

Subject to change without notice.
For additional information please contact
application@milestonesrl.com

General precautions

Always use hand, eye and body protection when operating with the microwave system.

